Lichtdusche

10006251

KLARSTEIN

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2014/30/EU (EMV)  
2014/35/EU (LVD)  
2011/65/EU (RoHS)  
93/42/EWG (Med)

Dear Customer,

Congratulations on purchasing this equipment. Please read this manual carefully and take care of the following hints to avoid damages. Any failure caused by ignoring the mentioned items and cautions mentioned in the instruction manual are not covered by our warranty and any liability.

## Technical Data

Item number	10006251
Power supply	220-240V, 50 Hz
Power consumption	2 x 36 W tubes
Light intensity	10.000 Lux

## Unpacking and checking


Every attempt is made to ensure your order is accurate and complete. However, to be sure that your order is correct, verify the contents of the box against your packing slip. The appliance is designed for simplicity of operation and set-up and requires minimal assembly. The following items are included your box:

- 1 (one) appliance
- 1 (one) Operator's Manual

Check the components according to the packing list. · Check for any damage or defects. Do not attempt to setup the appliance if anything is damaged or defective. Contact your supplier immediately if anything is found to be damaged or defective.

## Product Description

- 1 Screen
- 2 Brackets
- 3 Mounting hole
- 4 Carrying handle
- 5 Fold-out stand
- 6 On/off switch
- 7 Retaining clips
- 8 Reflector
- 9 Tubes, 2 x 36 W
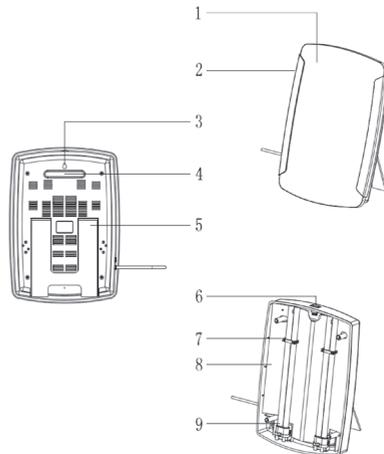
# About the Daylight

## Instructions for Use

Congratulations on the purchase of your new “Daylight”. The product will offer many years of reliable usage if handled correctly. Please keep your receipt, these instructions and the original packaging safe in case of future requirement. Please read the instructions carefully to gain the most benefit from your new “Daylight”, although the usage guidelines in this document do not replace the advice of your qualified medical advisor.

## Intended Use:

Medical intended use of this appliance is to treat SAD (seasonal affective disorder) or winter depression. “Daylight” is intended to boost people energy level, to adjust circadian rhythms and make people feel more awake. “Daylight” is safe and medically proven to relief from the winter blues.

## How does the “Daylight” work?

In autumn and winter, the seasons with the least sunlight because of shorter days, the lack of sunlight and a problem with certain brain chemicals stops the hypothalamus working properly. The lack of light is thought to affect the production of the hormone melatonin (that makes you sleepy). “Daylight” uses very bright lights to simulate the sunlight during the darker autumn and winter, whilst eliminating the dangerous UV rays contained in sunlight. The additional light helps your brain to reduce to production of hormone melatonin, adjust the levels of melatonin that can help the symptoms of SAD or winter depression.

## Seasonal Affective Disorder:

Seasonal Affective Disorder (SAD) as known as winter depression or winter blues. It may cause people:

- Feel low, listless and fatigued
- Need for more sleep
- Have trouble concentrating
- Less active than normal
- Lacking in energy

## Contraindication:

The following conditions may be contraindication for use of this appliance:

- suffer(ed) from severe depression
- recent eye surgery or a diagnosed eye condition which your doctor advised you to avoid bright light on
- specific medication (i.e. certain antidepressants, psychotropic drugs or malaria tablets);
- have high blood pressure
- skin sensitive to light, such as systemic lupus erythematosus

Always consult your doctor before your start using this appliance if you have any of the above conditions.

### Possible side effects:

The following side effects you may experience:

- headaches
- Eyestrain
- sleeping problems
- Nausea

### Measures to avoid, minimize or alleviate these side effects:

- Increase the treatment distance
- Avoid using the appliance during the late evening or before sleep
- Stop using the appliance for a few days to make the side effects disappear and then try again.

## Signs and Symbols

The following symbols appear in this manual, on the appliance:

 WARNUNG	Warning instruction indicating a risk of injury or damage to health		Read operating instructions before use
 WICHTIG	Safety note indicating possible damage to the unit/ accessory		Consult Operator's Manual
	Note on important information		Caution: Consult accompanying documents
O/I	OFF/ON		Class II Equipment
SN	Specifies serial number of the "Daylight"		Date of Manufacture
	Authorized Representative in the European Community		Manufacturer
IP20	Ingress Protection Rating		Keep Dry

## Safety and Performance Standards

This appliance has been designed and manufactured to conform to the highest standards of safety and performance including electromagnetic compatibility (EMC). This appliance conforms to the applicable requirements of the following:

This appliance meets the requirement of Medical Device Directive 93/42/EEC for medical appliances.

## Normative references

Guidance and manufacturer's declaration – electromagnetic immunity			
The appliance is intended for use in the electromagnetic environment specified below. The customer or the user of the appliance should ensure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment guidance
Electrostatic discharge (ESD) IEC 61000-4-2	6 kV contact 8 kV air	6 kV contact 8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrical fast transient/burst IEC 61000-4-4	2 kV for power supply lines 1 kV for input/output lines	2 kV for power supply lines 1 kV for input/output lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	1 kV differential mode 2 kV common mode	1 kV differential mode 2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	<5 % UT (>95 % dip in UT) for 0,5 cycle 40 % UT (60 % dip in UT) for 5 cycles 70 % UT (30 % dip in UT) for 25 cycles <5 % UT (>95 % dip in UT) for 5 sec	<5 % UT (>95 % dip in UT) for 0,5 cycle 40 % UT (60 % dip in UT) for 5 cycles 70 % UT (30 % dip in UT) for 25 cycles <5 % UT (>95 % dip in UT) for 5 sec	Mains power quality should be that of a typical commercial or hospital environment. If the user of the appliance requires continued operation during power mains interruptions, it is recommended that the appliance be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE: UT is the a.c. mains voltage prior to application of the test level.			

### Guidance and manufacturer's declaration – electromagnetic immunity

The appliance is intended for use in the electromagnetic environment specified below. The customer or the user of the appliance should ensure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment guidance
Conducted RF IEC 61000-4-6 Radiated RF IEC 61000-4-3	3 Vrms 150 kHz to 80 MHz 3 V/m 80 MHz to 2.5 GHz	3V  3 V/m	Portable and mobile RF communications equipment should be used no closer to any part of the appliance, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.  Recommended separation distance  $d=1.2\sqrt{P}$ $d=1.2\sqrt{P}$ 80MHz to 800MHz $d=2.3\sqrt{P}$ 800MHz to 2.5 GHz  where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in metres (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, a should be less than the compliance level in each frequency range.b Interference may occur in the vicinity of equipment marked with the following symbol:  

**NOTE 1:** At 80 MHz and 800 MHz, the higher frequency range applies.

**NOTE 2:** These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

a	Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the appliance is used exceeds the applicable RF compliance level above, the Medical Daylight should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the appliance.
b	Over the frequency range 150 kHz to 80 MHz, field strengths should be less than [V] V/m.

Guidance and manufacturer's declaration – electromagnetic emissions		
The appliance is intended for use in the electromagnetic environment specified below. The customer or the user of the appliance should ensure that it is used in such an environment.		
Emissions test	Compliance	Electromagnetic environment – guidance
RF emissions CISPR 11	Group 1	The appliance uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	The appliance is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Complies	

Recommended separation distances between portable and mobile RF communications equipment and the Medical Daylight			
The appliance is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the Medical Daylight can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the appliance as recommended below, according to the maximum output power of the communications equipment.			
	Separation distance according to frequency of transmitter m		
Rated maximum output power of transmitter W	150 kHz to 80 MHz	80 MHz to 800 MHz	800 MHz to 2,5 GHz
	$d = \left[\frac{3,5}{V_1}\right]\sqrt{P}$	$d = \left[\frac{3,5}{E_1}\right]\sqrt{P}$	$d = \left[\frac{7}{E_1}\right]\sqrt{P}$
0,01	0,12	0,12	0,23
0,1	0,38	0,38	0,73
1	1,2	1,2	2,3
10	3,8	3,8	7,3
100	12	12	23
For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.			
<b>NOTE 1:</b> At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.			
<b>NOTE 2:</b> These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

# Safety Instructions

- Before use, ensure that all packaging materials are removed and that there is no visible damage to the appliance. When in doubt, do not use the appliance and contact your dealer or the customer service address provided.
- The appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Keep the cord away from hot objects and open flames.
- The mandatory minimum distance between the appliance and flammable objects is 1.0 meter.
- Connect the appliance only to the mains voltage listed on the type plate.
- Always remove the mains plug and let the appliance cool off before you touch it.
- Do not touch the appliance with wet hands while it is plugged in; do not allow any water to be sprayed onto the appliance. The appliance must be operated only when it is completely dry.
- Protect the appliance from strong impacts.
- Avoid pulling the cord when removing the plug from the socket. Instead, hold the plug.
- If the cord or the housing is damaged, contact customer service or the dealer, since special tools are needed to make repairs. Never try to repair the appliance yourself.
- Disconnection from the mains supply is only guaranteed if the mains plug has been removed from the socket.
- The appliance is not suitable for use in the presence of flammable anesthetic mixtures with air or with oxygen or nitrous oxide.
- The appliance requires no calibration, preventive inspection and maintenance.
- The appliance is not repairable and contains no user serviceable parts.
- The operator shall not touch the appliance and the patient simultaneously.
- No modification of this appliance is allowed. Do not modify this equipment without authorization of the manufacturer.
- If this appliance is modified, appropriate inspection and testing must be conducted to ensure continued safe use of the equipment.

## Important Hints on Use

- This appliance is not intended for commercial or clinical use, but only for individual, private household use.
- If you are taking medication such as pain relievers, antihypertensives or antidepressants, consult your doctor before using the appliance.
- People with retinal diseases as well as diabetics need to be examined by an ophthalmologist before beginning to use the appliance. Do not use if you have an eye disease such as cataracts, glaucoma, diseases of the optic nerve in general and inflammation of the vitreous humour.
- If you have health concerns of any kind, contact your general practitioner!

## Repairs

- Only open the appliance to replace the tube. Please do not attempt to repair the appliance yourself. This could result in serious injuries. Failure to observe this regulation shall void the warranty.
- If this appliance is damaged in any way or it is not working, it should not be operated. Please try to fix the problem according to FAQ table. If the problem still not resolved, please contact customer service or the dealer, since special tools are needed to make repairs. Never try to repair the appliance yourself.

# Use and Operation

1. Set the daylight bulb on a stable surface, such as a table, and stabilize by using the folding tripod.
2. Plug the power cord into the socket and press the on-off switch at the top of the lamp.
3. After use, press the on-off switch at the top of the lamp, then unplug the power cord from the socket.
4. Connect the unit only to the mains voltage listed on the type plate.

## Fastening of Daylight on the wall

There is a mounting hole on the back side of the appliance. At this mounting hole the appliance can be fixed on the wall with screw. Depending on the type of wall, an appropriate screw must be used, may be with pegs. E.g. for concrete wall, M4, 25mm long Stainless Steel Self-Tapping Screw which is screwed into a anchor. The screw and the anchor are not included in the scope of delivery. If the mounting hole or back enclosure is cracked or damaged. Do not fix the appliance on the wall. Please contact customer service or the dealer for assistance.

## Getting Started

- Press the ON/OFF button of the appliance to switch on the light.
- Sit at a distance of 20 to 50cm (from the face) from the lamp. During this time, you are free to do everyday activities like reading, writing and making phone calls.
- Look directly into the light intermittently, for short periods of time.
- Do not look directly into the light for the entire period of use. This could overstimulate the retina.
- Start with a shorter exposure time, and gradually increase over a week.

**Note:** The first applications could cause eyestrain and headaches, however these will lessen over subsequent sessions as the nervous system adapts to the new stimuli. In this case abort the use and reduce the usage time thereafter. Do not touch the appliance in the duration. Let the appliance cool down for at least 10 minutes before reposition it.

## Usage times:

- The recommended distance is between 20 and 50cm from the lamp.
- At 20cm distance, the lamp can be used for up to 0,5 hour.
- At 40cm distance, the lamp can be used for up to 1 hour.
- At 60cm distance, the lamp can be used for up to 2 hours.
- Basically - the closer you are to the lamp, the shorter the duration of use will be.

## Length of use

Repeat use during low-light seasons/low light conditions for at least 7 days consecutively. You made require longer than 7 days, depending on individual needs. Ideally, the lamp should be used in the morning.

## Cleaning and Care

Before cleaning, the lamp must be switched off and given time to cool down. To clean, use a damp cloth and some detergent if necessary. Pay particular attention that no water penetrates into the interior of the lamp. Do not use harsh cleaning agents and never put the lamp into water. The lamp cannot be cleaned in a dishwasher. Once connected, do not touch the lamp with wet hands or expose to water in any way. The lamp must be dry when used.




## Replacing the Lamp Tubes

In case the lights are faulty, or when they come to end of their life, the two light tubes can be replaced. NOTE: Do not replace the lamp tubes which are not provided by the supplier. To do so, proceed as follows:

### Important:

Always unplug the appliance and let it cool off.

For safety reasons, do not remove any components or screws other than those specified below.
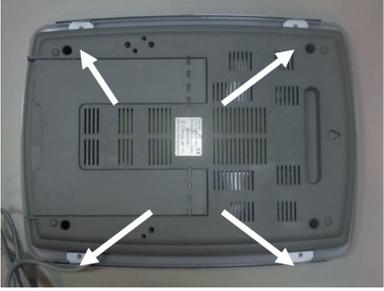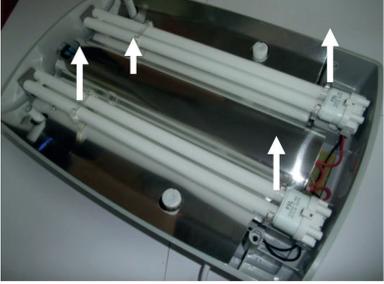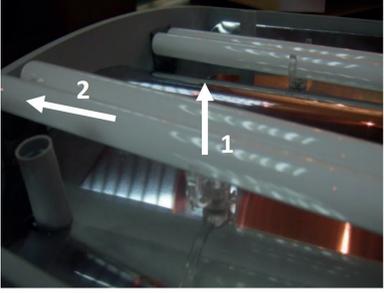
<p>Slacken the retaining screws on the back using a screwdriver, lay the apparatus on its back and remove the screen.</p> <p><b>Note:</b> The two brackets on the front of the screen will become loose when the screws are slackened. Take care not to lose them, or the screen cannot be refitted correctly.</p>	 A top-down view of the back of the lamp housing. Four white arrows point to screws located at the corners of the back panel, indicating they should be loosened.
<p>Cut the cable ties and bend the tube's four transparent retaining clips by pulling them gently upwards.</p> <p><b>Note:</b> Cable ties are used for transport protection only.</p>	 An interior view of the lamp housing showing two fluorescent tubes. Four white arrows point upwards to the transparent retaining clips that hold the tubes in place.
<p>Carefully pull the tube up over the housing rim (1) and out of the housing, then sideways out of the frame.</p>	 A close-up view of a fluorescent tube being pulled out of the housing. A white arrow labeled '1' points upwards to the rim of the housing, and another white arrow labeled '2' points to the left, indicating the direction to pull the tube out.
<p>When inserting the tube, following these instructions in reverse order.</p>	

## Cleaning Up a Broken Tube

The light tubes contain a small amount of mercury sealed within the glass tubing, when a tube breaks in your home, some of this mercury is released as mercury vapor. To minimize exposure to mercury vapor, we recommend that users follow the cleanup and disposal steps as follows:

### Before Cleanup:

Have people and pets leave the room and air out the room for 5-10 minutes by opening a window or door to the outdoor environment. Shut off the central forced air heating/ air-conditioning system, if you have one. Collect material needed to clean up broken tube, just like: stiff paper or cardboard, sticky tape, damp paper towels or disposable wet wipes.

### During Cleanup:

Vacuuming is not recommended unless broken glass remains after all other cleanup steps have been taken. Vacuuming could spread mercury-containing powder or mercury vapor. Be thorough in collecting broken glass and visible powder. Place cleanup materials in a sealable container.

### After Cleanup:

Promptly place all bulb debris and cleanup materials outdoors in a trash container or protected area until materials can be disposed of properly. Avoid leaving any bulb fragments or cleanup materials indoors. If practical, continue to air out the room where the bulb was broken and leave the heating/air conditioning system shut off for several hours.

## Storage

When the appliance is not being used, it should be stored in a dry location out of reach of children. Do not keep the it connected to a power source. The unit is designed to be stored in the following environmental conditions:

- Ambient temperature from 5-50 °C.
- Relative humidity of max. 90%.

## Troubleshooting

Trouble	Possible Cause	Solution
No light.	Main Power OFF.	Turn on the switch.
	No Power.	The appliance is not correctly connected to your power supply is defective.
	Loose Light tube.	Re-install the light tube.
One side of the product is out of order.	Defective light tube.	Change the light tube.
	Light tube at end of life.	Change the light tube.
	Loose Light tube.	Re-install the light tube.

This chapter lists the questions most frequently asked about the appliance. If you cannot find the answer to your question, please contact customer service or the dealer for assistance.

## Technical Data

Item number	10006251
Dimensions (W x H x D)	340 x 480 x 80 mm
Weight	2 kg
Lighting elements	2 x 36 W tubes (2G 11), Daylight (full spectrum)
Nominal data for replacement lighting elements (provided by light tube supplier)	FSL - YDW36-H.RR 36W/865 2G11 4pin / Philips - MASTER PL-L 36W/865/4P Service Life: 8000 hours
Fuse	T3.15AH 250V
Power	72 W
Light intensity	> 10.000 lux in 15 cm
Mains connection	220-240 V - 50 Hz
Operating Environment	50 - 77 °F/ 10 - 25 °C, and relative humidity from 30% to 90%. Pressure altitude: <2000 m Atmospheric pressure: 700-1060 hPa
Storage Environment	41 - 122 °F/ 5 - 50 °C, and relative humidity of max. 90%. Pressure altitude: <2000 m Atmospheric pressure: 700-1060 hPa
Transport Environment	41 - 122 °F/ 5 - 50 °C, and relative humidity of max. 90%. Pressure altitude: <2000 m Atmospheric pressure: 700-1060 hPa
Product Classification	1. Class II equipment 2. No applied part 3. Protection against ingress of water: IP20 4. Not category AP / APG equipment 5. Mode of operation: Continuous operation.

Subject to technical modification.

## Hints on Disposal

According to the European waste regulation 2002/96/EC this symbol on the product or on its packaging indicates that this product may not be treated as household waste. Instead it should be taken to the appropriate collection point for the recycling of electrical and electronic equipment. By ensuring this product is disposed of correctly, you will help prevent potential negative consequences for the environment and human health, which could otherwise be caused by inappropriate waste handling of this product. For more detailed information about recycling of this product, please contact your local council or your household waste disposal service.

## Declaration of Conformity

European Authorised Representative:

I-Plus Company Ltd.  
4B Christchurch House, Beaufort Court,  
Medway City Estate, Rochester, Kent. ME2 4FZ, UK  
Tel: (44)-20-71936086  
Email: [i-plus@zeroplus.hk](mailto:i-plus@zeroplus.hk)

**Manufactured for:** Chal-Tec GmbH, Wallstraße 16, 10179 Berlin, Germany.

This product is conform to the following European Directives:

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**Manufactured for:** Chal-Tec GmbH, Wallstraße 16, 10179, Berlín (Alemania).